

5" touchscreen
individual display
the maximum comfort of operation



www.radwag.com

X2

The new series
of laboratory balances

ADVANCED LEVEL

X2

The new series
of laboratory balances
RADWAG
2015

Comfort of operation due to 5" capacitive touchscreen

The X2 series balances feature the latest generation capacitive display providing the maximum comfort of use, available right at your fingertips. Ease of operation, clear menu and practical arrangement of the display guarantee the best ergonomics for your everyday tasks.

Personalization taken to the next level

The X2 is the only balance available on the market that, provides a user with the option to design display using wide selection of widgets. Customised display of X2 series offers direct access to the results of your work and other important information directly from the home screen.

Unlimited possibilities of communication

A wide array of available interfaces facilitate selection of the most optimal means for communication. The X2 series balances offer unlimited possibilities for cooperation with external devices, providing printing, copying, archiving and data transfer.

Ergonomics of operation

Built-in IR sensors allow numerous operations (e.g. tarring, transmitting the result to a printer or selecting successive steps of a particular process, etc.) to be performed hands-free, by simply moving a hand across the sensor.

Data safety


The X2 series balances feature ALIBI memory that is a warranty for safety and automatic recording of your measurements. Options such as data preview, copying and archiving are available to users, as well.

Weighing data management

USB interface facilitates quick transfer and copying of any results of your work (measurements, reports, databases) to other balance.



AS X2
Analytical balances

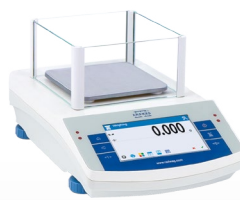

PS X2
Precision balances



PS 10 X2
Precision balances



APP X2
Precision balances

Max capacity [Max]	52 - 310 g	200 - 6000 g	1 - 10 kg	10 - 32 kg
Readability [d]	0,01 - 0,1 mg	1 - 10 mg	0,01 - 0,1 g	0,01 - 0,1 g
Display	5" colour capacitive touchscreen	5" colour capacitive touchscreen	5" colour capacitive touchscreen	5" colour capacitive touchscreen
Adjustment	internal automatic	internal automatic	internal automatic	internal automatic
Stabilization time	3,5 - 6 s	1,5 - 2 s	3 s	2 - 3 s
Pan size	ø 70 mm or ø 85 mm	128×128 or 195×195 mm	195×195 mm	348×260 mm
Interfaces	USB A, USB B RS 232 (2 ports) Wi-Fi (option)	USB A, USB B RS 232 (2 ports) Wi-Fi (option)	USB A, USB B RS 232 (2 ports) Wi-Fi (option)	USB A, USB B RS 232 (2 ports) Wi-Fi (option)

